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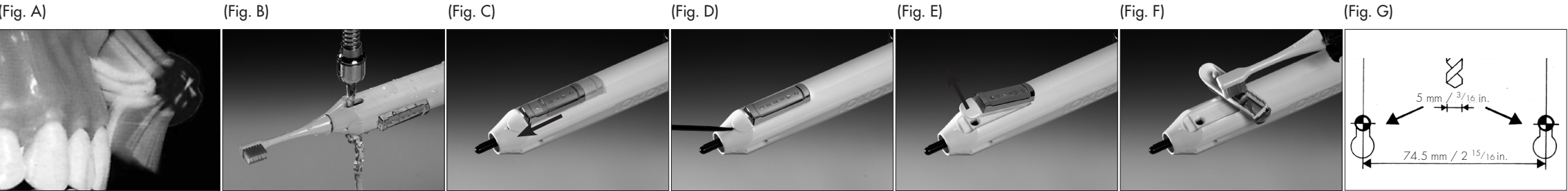
Using your BROXO® OraBrush™

- Connect the step-down converter to the appropriate outlet with the switch on position O.
- Completely insert a stem brush until an audible “click” is heard
- Run water on the stem brush bristles before applying toothpaste, then insert the stem brush in the mouth
- Start the device by sliding the switch to the desired power 1,2 or 3.
- After each use, rinse the device well under running water (Fig. B).
- It is suggested to hand the device back on the wall holder after use so that water may drain and it may dry optimally. If the wall is not used, the device may be set on its base and the stem brushes hung from the adhesive mini wall holder (not included, available on order).
- It is recommended to pull the step-down converter from the outlet after use.

You will quickly realize that you now own the best performing prophylaxis device for teeth and gums. Here are some instructions on how to keep it in perfect working order and insure continued best results.

Caring for your BROXO® OraBrush™

The entire device can be cleaned under lukewarm water complemented with soap or a non-damaging cleaner compatible with plastic materials. In case of calcium deposits, only clean the device with water mixed with white vinegar.

(Fig. B) So as to prevent dirt buildup (toothpaste, calcium or other), it is essential to properly rinse the device under lukewarm running water, and particularly to let water flow through the 2 elongated openings on each side of the device, as well as on the switch. Feel free to let the device run and to activate the switch while rinsing for a more thorough clean. If toothpaste seeps into the switch, it is easy to take it apart for cleaning:

- (Fig. C) Slide it forward to position 3
- (Fig. D) Insert a pointed, non-sharp object (e.g. knitting needle, plastic stylus) in the notch at the tip of the switch and use it as a lever to pop it out.
- (Fig. E+F) Remove the front of the switch. Clean the front of the device and the switch with an old stem brush or equivalent. Make sure not to damage the internal black rubber protection that can be seen through the 2 elongated openings.
- To replace the switch back into the device, set it to position 3 and insert the plastic pin at the back of the switch into the notch at the base of the switch in the device casing. Press on the front of the switch, an audible “click” will signal its proper attachment.

For an optimal effectiveness, it is recommended to replace the stem brushes every 3 months. For a complete oral hygiene it is recommended to complement the use of the BROXO® OraBrush™ with a BROXO® OraJets™ irrigator.

Installing the wall holder

- Pick the desired spot
- Drill 2 holes as shown in Fig. G.
- Insert the plastic screw holders and drive in the screws, leaving enough of their heads out to slide the wall holder onto them.

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Place the stem brush between teeth and gums (fig. A)

Type of use	Age categories	Frequency	Angle adjustment position									
			1 st day	4 th day	7 th day	1 th day	15 th day	21 th day	30 th day	45 th day	+60 th day	
1. Normal usage	Children	Twice daily	1			2						
	Adults	1-3 minutes	2			3						
2. Gingivitis Post Gingivectomy		3x3 min 3x2 min	1			2			2		3	
3.Parodontitis		3x4-5 min 3x3 min 3x2 min 2x2-3 min	1			2						
								2		3		
then normal usage												
4. Wearers of orthodontic appliances	Children	Twice daily	1				2					
	Teenagers Adults	2-3 minutes	2					3				

In special cases it is recommended the dentist prescribes the utilisation procedure

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Safety instructions

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- The external flexible cable or cord of this transformer cannot be replaced; if the cord is damaged, the transformer shall be scrapped.

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